


Carebidet Quick Guide 12

Cleaning & Routine Maintenance

Disinfecting the Water Tank and Reservoir

Use this with the training video

This quick guide summarizes the procedure shown in the approved Carebidet training video. For complete operating instructions and cautions, use the Carebidet User Manual and the full video.

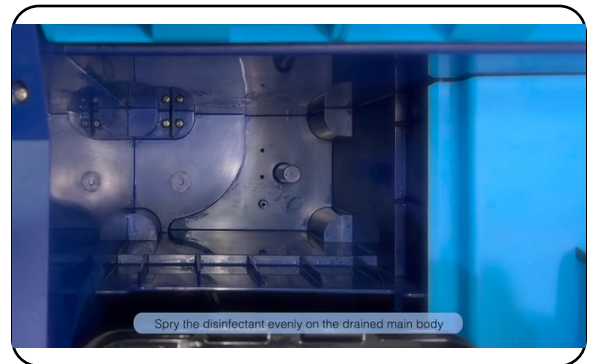

Quick steps

1. Open the main body and remove the water tank and water filter as shown.
2. Connect the drain hose to the water outlet in the main body.
3. On the control panel, select Cleaning Settings > 1-015 Manual Drainage and press OK to drain the system.
4. After drainage, apply the approved disinfectant evenly to the drained main-body reservoir surfaces.
5. Wipe the main body thoroughly with a dry cloth after the disinfecting step.
6. Reinstall the water filter in the correct direction, placing the center hole over the water-filter guide first.
7. Reinstall the water tank and close the main-body cover.

Important

Stop and re-check the setup if the individual has pain, skin pinching, leakage, an unresolved alert, or an insecure connection.